

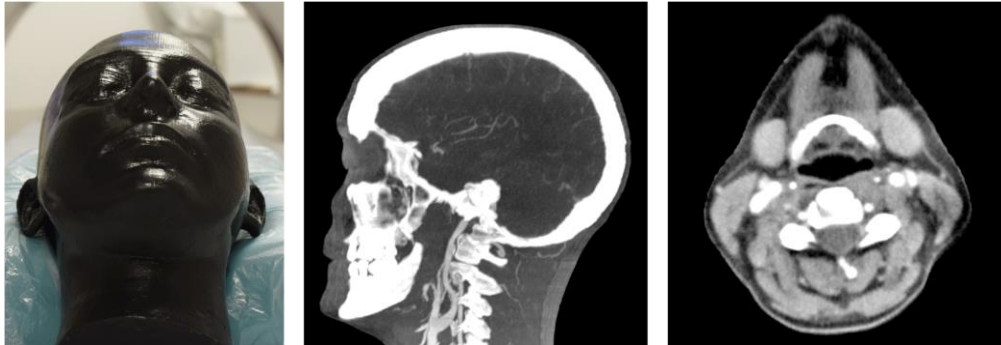


Head and Neck Phantoms

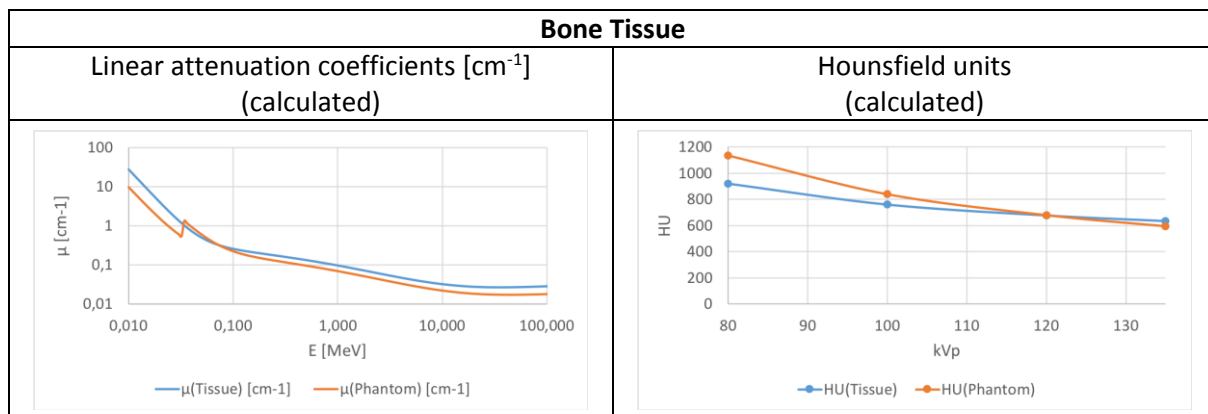

Size	Approx. 17 x 20 x 28 cm
Weight	Approx. 4 kg

Phantom attenuation mainly influenced by iodine. Attenuation values are calibrated to be equivalent to human tissue at 120 kVp. Calibration to other spectrum energies is possible upon request.

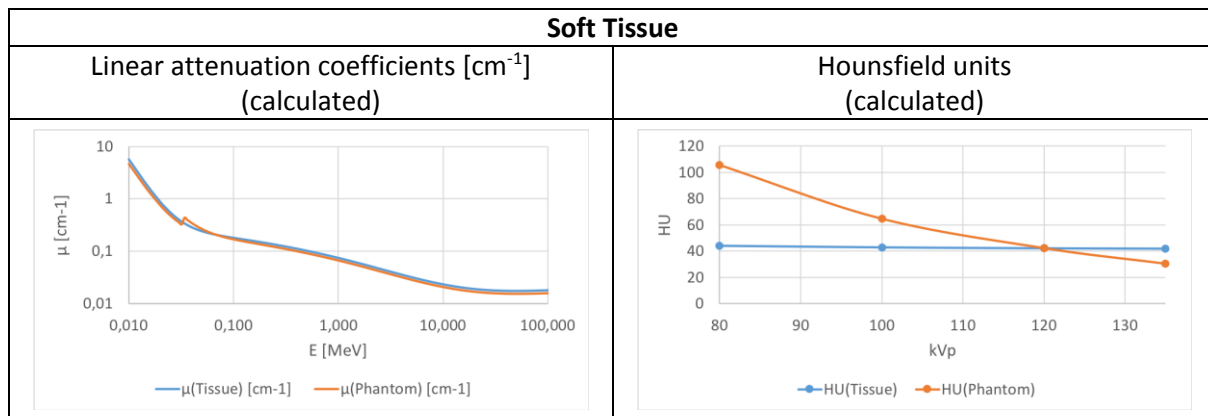

General indications

- Phantoms are manufactured of a cellulose-polymer composite material with similar properties to hardwood. If treated carefully, they will last for a long period.
- The phantoms are coated with a protective layer. If the protective layer is unharmed, the phantoms can be cleaned using a damp cloth (water or mild detergent).
- Protect from direct sunlight.
- Maintain a storage temperature of 10 °C to 30 °C.
- If the phantom is exposed to temperatures below 0 °C or above 60 °C, even for a short period of time, it can be seriously damaged.